

Murisphere Mobile Scan Tutorial

Professional Step-by-Step Guide for Printed Cage Cards

Murisphere Operations

2026-03-04

Contents

Purpose	1
Who Should Use This	2
Before You Start	2
Visual Walkthrough	2
Figure A - Printed Card to Phone Flow (Illustration)	2
Figure B - Troubleshooting Matrix (Illustration)	2
Step-by-Step: Phone Scan of Printed Card	2
Step 1 - Print Fresh Cards	2
Step 2 - Scan with Phone Camera	3
Step 3 - Validate Identity	3
Step 4 - Enter Edit Mode	3
Step 5 - Perform Update	3
Required Quality Checks	3
High-Value Troubleshooting	3
If camera does not detect code	3
If code opens but wrong destination	3
If QR appears broken	4
Trainer Checklist	4
Daily Room SOP (Quick)	4

Purpose

This tutorial shows technicians exactly how to use a **phone camera** to scan a **printed cage card** and open the correct cage record in Murisphere.

Who Should Use This

- Technicians working in-room
- PI delegates verifying cage records
- Facility trainers onboarding new users

Before You Start

1. Confirm Murisphere is reachable from your phone browser.
2. In the Cages tab, set Scan Base URL to a phone-reachable address.
3. Print cards using Generate + Print.

Visual Walkthrough

Figure A - Printed Card to Phone Flow (Illustration)

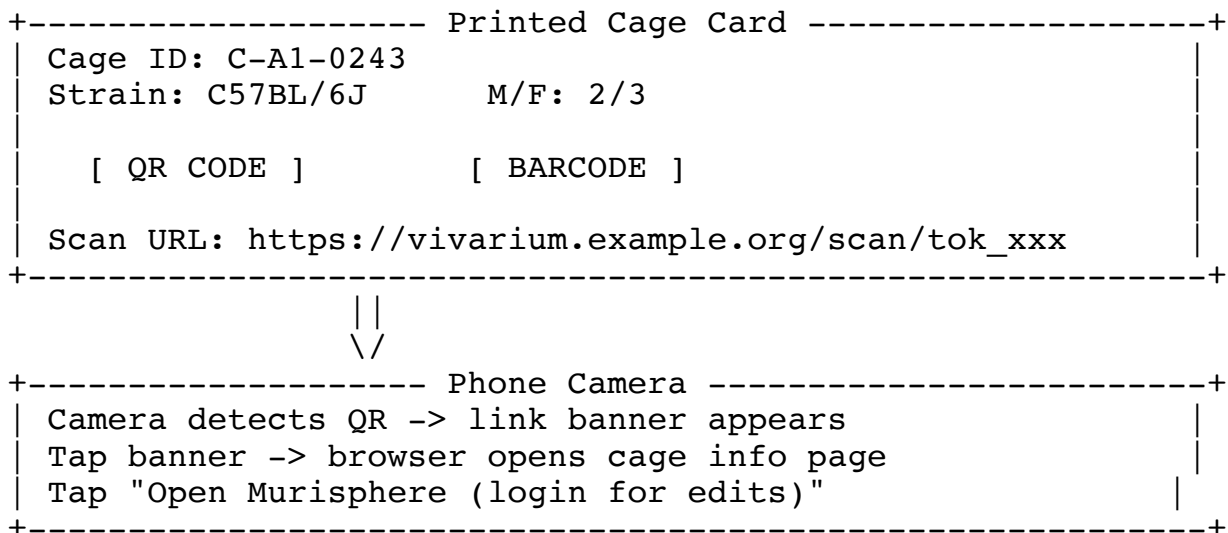


Figure B - Troubleshooting Matrix (Illustration)

Issue	Likely Cause	Action
No scan prompt	Low contrast / glare	Improve lighting
Wrong page opens	Bad Scan Base URL	Use LAN IP/domain
Blank QR or broken image	Cached JS / blocked lib	Hard refresh, re
Slow open	Weak room network	Retest on strong

Step-by-Step: Phone Scan of Printed Card

Step 1 - Print Fresh Cards

1. Open Cages.
2. Confirm Scan Base URL points to a phone-reachable host.
3. Click Generate + Print.

4. Verify each card has a visible square QR.

Step 2 - Scan with Phone Camera

1. Hold the printed card under even lighting.
2. Open default camera app:
 - iPhone: Camera app
 - Android: Camera app
3. Point camera at the **QR code** until a link banner appears.
4. Tap the banner/link.

Step 3 - Validate Identity

1. Confirm cage ID on the opened page matches printed cage ID.
2. Confirm location and M/F counts are plausible.
3. If mismatch: stop and rescan before editing.

Step 4 - Enter Edit Mode

1. Tap Open Murisphere (login for edits).
2. Sign in if prompted.
3. App should land in scan/edit context for that cage.

Step 5 - Perform Update

1. Update counts/status/notes.
2. Save changes.
3. Confirm save acknowledgement and audit trail entry.

Required Quality Checks

- Do not update if cage ID mismatch.
- Do not use localhost for printed scan URLs.
- Reprint cards after base URL changes.

High-Value Troubleshooting

If camera does not detect code

- Move 10-20 cm from card.
- Avoid reflections and shadows.
- Increase print quality.

If code opens but wrong destination

- Check Scan Base URL in Cages tab.
- Use LAN IP or production DNS.

If QR appears broken

- Hard refresh page (Cmd+Shift+R / Ctrl+F5).
- Regenerate cards.
- If still broken, capture browser console errors and report with timestamp + cage ID.

Trainer Checklist

1. Demonstrate one complete scan -> edit -> save cycle.
2. Have trainee repeat with a second cage.
3. Verify trainee can troubleshoot URL and scan issues.
4. Confirm trainee understands ID verification before edits.

Daily Room SOP (Quick)

1. Check scan URL base setting.
2. Print cards for active cages.
3. Run scan-update rounds.
4. Review scheduled tasks and compliance alerts.