

MURISPHERE TRAINING

Mobile Scan Tutorial (Printed Cage Card)

Professional phone-first workflow with vectorized screenshots for crisp print and PDF output.

Goal

Use a phone camera to scan a printed cage card and open the correct Murisphere cage record safely and quickly.

Screenshot 1 - Secure Login

Murisphere Secure Login

Tip
Use technician account for in-room updates.

Sign In
Role-based access with full audit trail

Email
tech@murisphere.local


Password

Sign In

Vectorized screenshot: Login panel on mobile-target web UI


Sign in with technician credentials before running room updates.

Screenshot 2 - Generate and Print Cage Cards

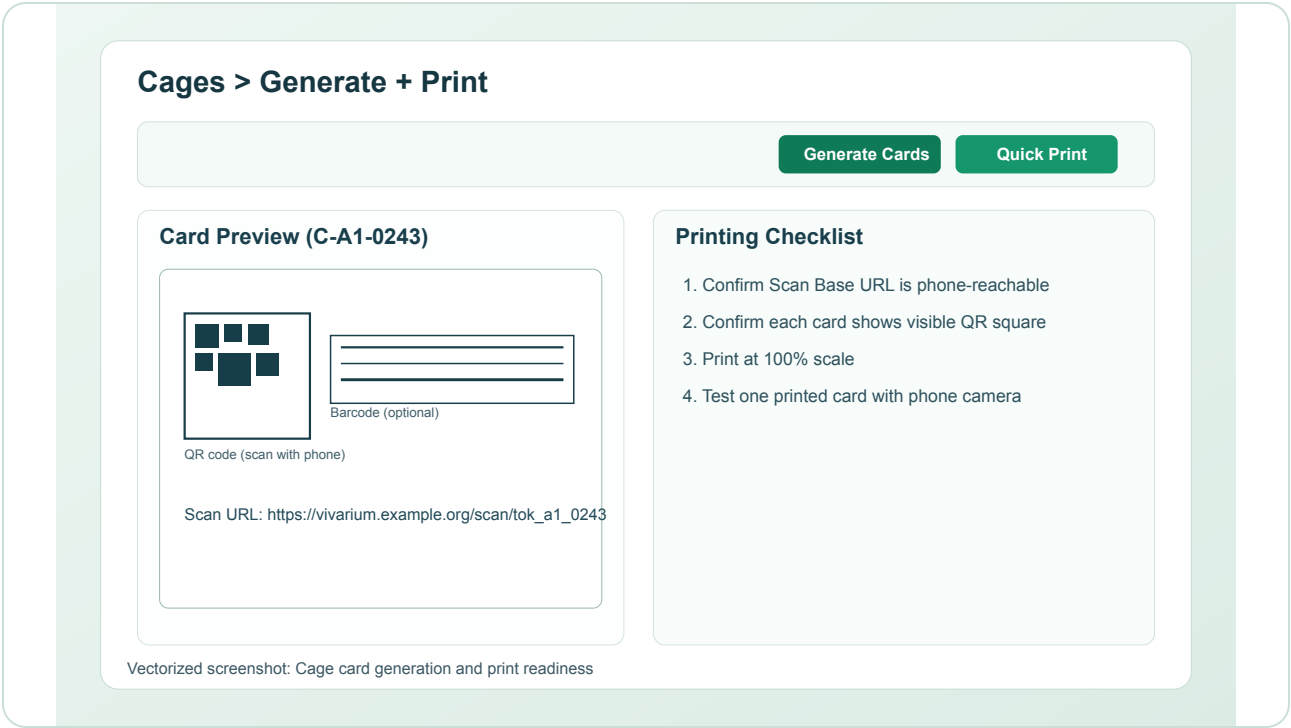

Confirm each card shows a visible QR block and print at 100% scale.

Screenshot 3 - Scan/Edit Workflow

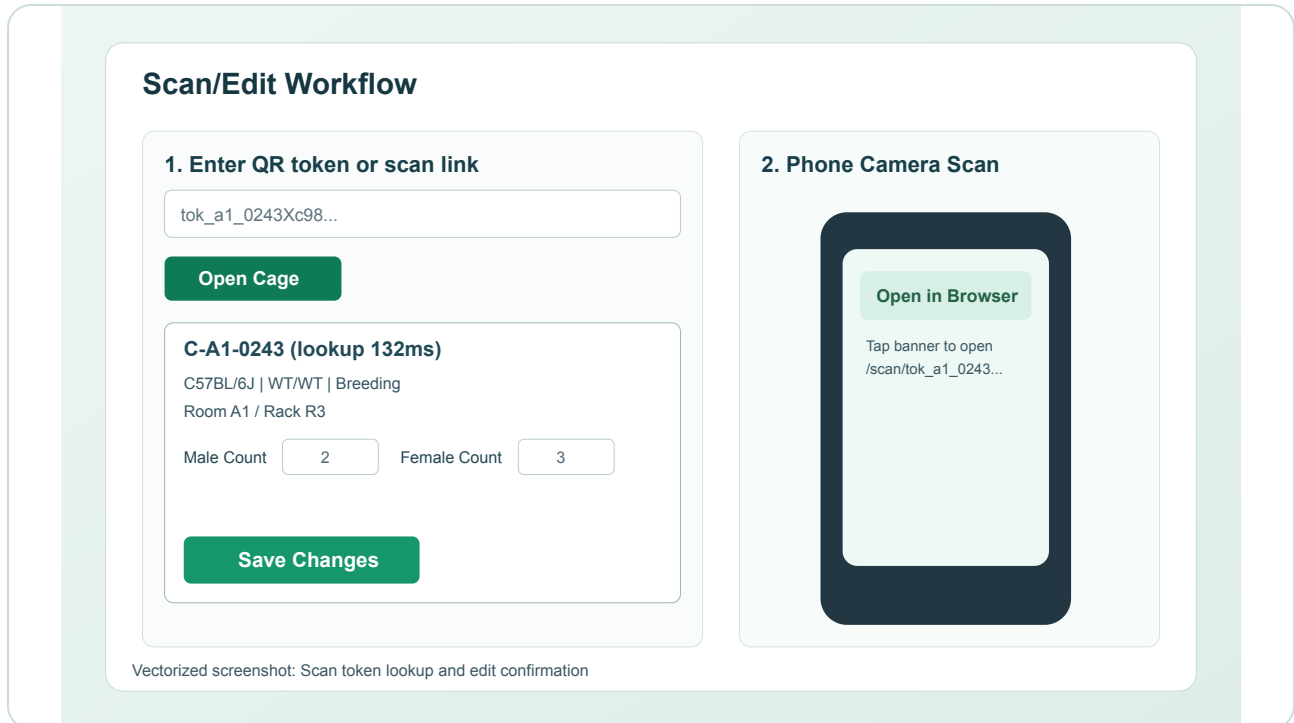

Phone camera opens the scan URL, then technician validates cage ID before saving edits.

Step-by-Step: Phone Scan of Printed Card

1. In **Cages**, set **Scan Base URL** to a phone-reachable host (LAN IP or production domain).
2. Click **Generate + Print** and verify the printed card contains a solid QR square.
3. Open the phone camera and point it at the printed QR code.
4. Tap the detected link banner to open cage info in browser.
5. Tap **Open Murisphere (login for edits)**, verify cage ID, then save updates.
6. If Wi-Fi drops, save actions are queued and auto-retried on reconnect. Confirm retry sync before ending rounds.

Critical Safety Checks

- **ID first:** never edit unless scanned cage ID matches printed card.
- **Protocol hard stop:** if protocol is expired, editing is blocked until protocol assignment is corrected.
- **Queue awareness:** if you see a queued-save message, reconnect and verify data synced before leaving the room.

Troubleshooting

- **No scan prompt:** improve lighting and flatten the printed card.
- **Wrong page opens:** correct `Scan Base URL` (never use localhost for printed cards).
- **Broken/blank QR:** verify server endpoint `/api/assets/qrcode.png?v=test` loads.
- **Save blocked by protocol:** check protocol expiry in Compliance tab and reassign valid protocol.
- **Queued save message:** verify network, wait for auto-retry, then reload cage to confirm persisted values.
- **Slow open:** retest on stronger room Wi-Fi.

Trainer Sign-Off

1. Technician demonstrates full scan -> verify -> edit -> save flow on two cages.
2. Technician demonstrates handling of one queued-save scenario (network drop simulation).
3. Technician identifies and escalates a protocol-expired block correctly.